

Precision Zero Air 18L


Engineered specifically to provide for flame support

Part number: 60-0019 (110V) | 60-0018 (230V)



Description

The Precision Zero Air 18L is engineered specifically to provide flame support gas to multiple GC instruments. It provides up to 18 liters per minute of hydrocarbon free air, which is suitable to support your FID, FPD and NPD. Expanding Peak Scientific's Precision range, the Precision Zero Air 18L, has the same stackable form-factor and compact footprint found in the rest of the Precision Series, ensuring it integrates seamlessly with other Precision models to provide a convenient, reliable and safe gas solution for specific GC gas requirements.

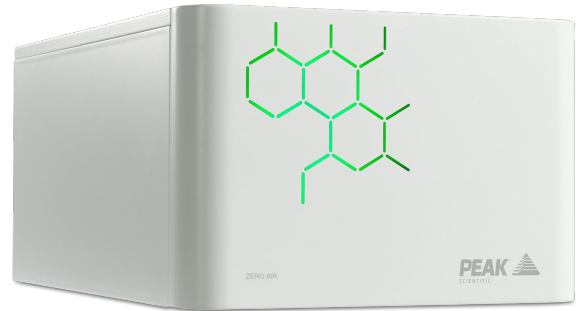

Application Types


GC Detector Gas

Generator Details

- **Gas types:** Air
- **Integrated compressor:** No

Key Features

- Integrates with existing Precision products, with same modular, stackable design
- Hydrocarbon content (as Methane) < 0.05ppm for market-leading purity
- Up to 18L/min of Zero Air Up to 80psi/5.5bar outlet pressure
- Minimum maintenance requirements, no expensive catalyst chamber replacements
- Greatly reduced risk of contaminants entering system and no risk of running out of gas
- Stylish and contemporary aesthetics and functional design
- Status indicative lighting feature (amber warmup, green ready)
- Compatible with Precision Air Compressor
- 12 month comprehensive on-site warranty

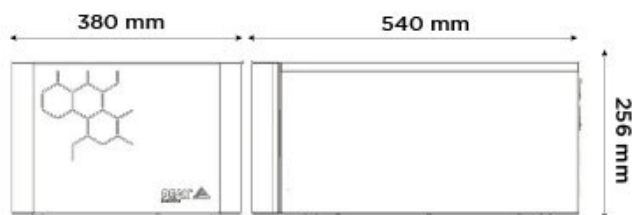

Technical Specifications

Precision Zero Air 18L Gas Generator

Gas Type	Air
Max Gas Flow:	18L/min
Max Output Pressure:	80psi/5.5bar
Hydrocarbon Content:	<0.05ppm
Max Output Pressure:	80psi/5.5bar
Gas Outlets Fitting:	1 x 1/4" BSPP
Min Inlet Pressure:	110psi/7.6bar
Max Inlet Pressure:	125psi/8.6bar
Min Air Inlet Flow:	18L
Min Inlet Air Quality:	ISO 8753-1:2010 Class 1.4.1 or Precision Air Compressor
Start Up Time:	60 mins
Power Consumption:	594 / 759 W
Voltage:	110 / 230 V
Frequency:	50 / 60 Hz
Current:	3.3 - 5.4 A
Heat output:	2590
Accreditations:	CE, FCC, CSA
Size (HxWxD) mm	256 x 380 x 540 mm
Size (HxWxD) Inches	10.1 x 15 x 21.3 inches
Generator Weight	25kg / 55lbs

[PEAK Protected]TM

Peak Scientific gas generators define the benchmark in reliability, convenience and performance in laboratories around the world, and come backed with a 12 month warranty. Beyond this period however you can ensure that your investment continues to be **[Protected]** by our comprehensive generator care cover.

Our world-class aftercare support packages deliver a program of scheduled preventative maintenance whilst giving you the reassurance of instant access to worldwide technical support and priority on-site response in the untimely event of a breakdown.

- For ordering parts visit: www.peakscientific.com/ordering/
- For service plans visit: www.peakscientific.com/service/service-plans/

Peak Scientific UK
Tel: +44 (0)141 812 8100
Fax: +44 (0)141 812 8200

Peak Scientific North America
Tel: +1-800-767-6532
Fax: +1-978-608-9503

Peak Scientific China
Tel: +86 21 5079 1190
Fax: +86 21 5079 1191
For a full list of our worldwide office locations, please visit:
Web: www.peakscientific.com
Email: discover@peakscientific.com

Peak Scientific's Quality Management System conforms to: ISO:9001

